
KYARA®

LTK600

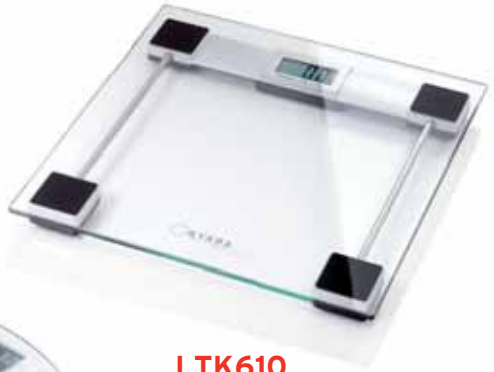LTK610

LTK615

BILANCIA PESAPERSONE

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MORETTI S.P.A.
Via Bruxelles 3 - Melegnano
52022 Cavriglia (Arezzo)
Tel. +39 055 96 21 11

Fax +39 055 96 21 200
www.morettispa.com
info@morettispa.com

INDICE

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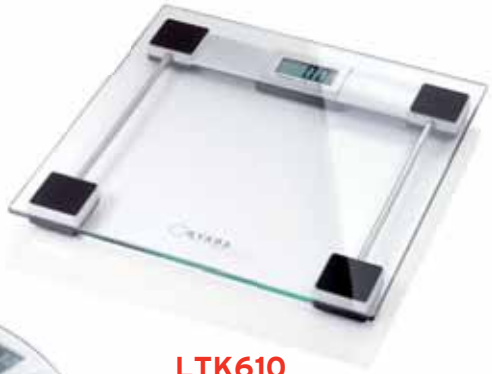
LTK610

LTK615

BATHROOM SCALE INSTRUCTION MANUAL

MORETTI S.P.A.
Via Bruxelles 3 - Melegnano
52022 Cavriglia (Arezzo)
Tel. +39 055 96 21 11

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CODES

- LTK600** Bathroom Scale with plastic platform
LTK610 Bathroom scale with square glass platform
LTK615 Bathroom scale with round glass platform

INTRODUCTION

Thank you for purchasing a **KYARA** bathroom scale by MORETTI S.p.A. This instruction manual contains advices and for your safety and for a correct use of the aid. We recommend reading thoroughly the instruction manual before using the bath lift. In case of doubts, contact your retailer. In case of doubts, we recommend to contact the dealer, who can help and advice you properly.

NOTE CHECK ALL PARTS TO VERIFY IT WON'T BE ANY SHIPPING DAMAGE. IN CASE OF SHIPPING DAMAGES, DO NOT USE. CONTACT THE DEALER FOR FURTHER INSTRUCTIONS. SHIPPING. IN CASE OF DAMAGE, PLEASE DO NOT USE THE DEVICE AND CONTACT YOUR DEALER FOR FURTHER INSTRUCTIONS.

INTENDED USE


This scale is intended to detect the weight of all the people, with only the maximum capacity limitation (see technical features).

GENERAL WARNINGS


- Do not use the product for a purpose not indicated in this manual. MORETTI S.p.A declines all responsibilities for any consequences resulting from an incorrect use of this product and from unauthorized alteration to the frame of the product. The manufacturer can modify the device or the instruction manual without prior notification.
- Do NOT install or use the device before reading and understanding this manual. If you are unable to understand the warnings, cautions or instructions, contact a healthcare professional dealer or technical personnel before attempting to install this equipment to avoid injure or damage
- Do not step on the scales with wet or dump feet.

DESCRIPTION OF PARTS

LTK600

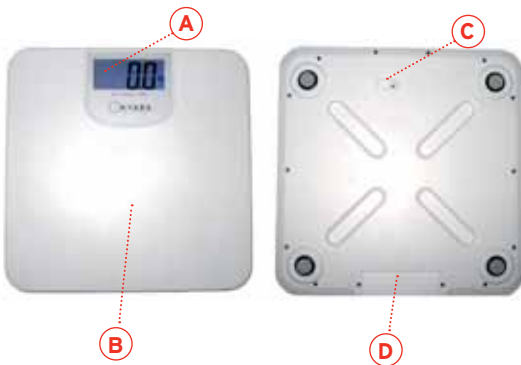- A** -Digital display with backlight
B -Footboard
C -Selector (Lb-Kg-St) =
 Pound-Kilogram-Stone
D -Battery compartment
 (4 batteries 1,5 V
 type AAA not included)

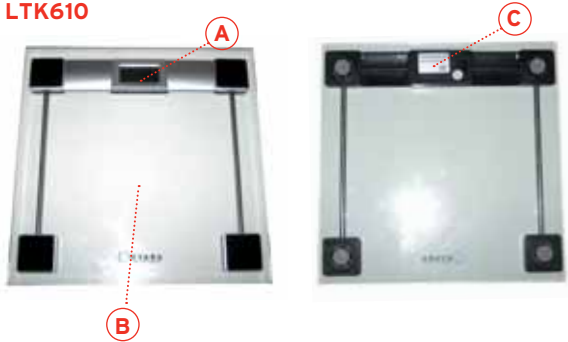LTK610

- A** - Digital display
- B** - Footboard
- C** - Battery compartment
(1 lithium battery 3V
CR2032 included)

NOTE INSIDE THE BATTERY COMPARTMENT, THERE IS A SELECTOR TO CHOOSE THE UNITS (LB-ST-KG) POUND-STONE- KILOGRAM.

LTK615

- A** - Digital display
- B** - Footboard
- C** - Battery compartment
(1 lithium battery 3V
CR2032 included)
- D** - Button to select the units
(Lb-St-Kg)
Pound-Stone Kilogram

ASSEMBLY INSTRUCTIONS

Battery installation

These kind of scales do not require special assembly operations, excluding the installation of the batteries that must be completed in the type and model stated in this manual and the polarity marked on the battery compartment.

HOW TO USE

LTK600

1. Install 4 "AAA" batteries in the battery compartment, matching the indicated polarity. Select the correct units on the back of the balance acting on the switch (LB = LB / KG = KG / ST = Stone) Please note that this must be done only once, during the first usage
2. Place your scale on a flat and clean surface
3. Step on the scale and put your feet evenly on the scales without moving
4. The back light will switch on and your correct weight will be displayed
5. Once it measured the weight, the display will flash three times and then it will continuously indicate the measured value

6. The weight will be displayed for about 5 seconds, then the scales will turn off in order to conserve battery power
7. To make an additional measurement, you have to get back again on the scale platform

LTK610

1. Install no.1 3V CR2032 lithium battery in the battery compartment, matching the indicated polarity
2. Select the correct units on the back of the balance acting on the switch that is always located inside the battery compartment (LB = LB / KG = KG / ST = Stone)
3. Please note that this must be done only once, during the first usage.
4. Place your scale on a flat and clean surface.
5. Switch on the scales by a slight tap on the glass platform.
6. When the display will show "0.0", step on the scale and put your feet evenly on the scales without moving
7. The weight will be displayed for about 5 seconds, then the scales will turn off in order to conserve battery power
8. To make an additional measurement, please repeat the procedure from point "4"

LTK615

1. Install no.1 3V CR2032 lithium battery in the battery compartment, matching the indicated polarity
2. Switch on the scales by a slight tap on the glass platform (LB= libbra / KG= Kilogrammo / ST= Stone)
3. Please note that this must be done only once, during the first usage
4. Place your scale on a flat and clean surface
5. When the display will show "0.0", step on the scale and put your feet evenly on the scales without moving
6. The display will switch on and your correct weight will be displayed.
7. Once it measured the weight, the display will flash three times and then it will continuously indicate the measured value
8. The weight will be displayed for about 5 seconds, then the scales will turn off in order to conserve battery power
9. To make an additional measurement, please repeat the procedure from point "4"

CLEANING

For the scale cleaning, use only a damp cloth and mild detergents. Then dry thoroughly before using the device.

CONDITIONS OF DISPOSAL

INSTRUCTIONS FOR THE PROPER DISPOSAL OF THE PRODUCT IN ACCORDANCE WITH DIRECTIVE 2002/96/EC

At the end of its life, the product should not be disposed of as waste. It can be taken to the special collection points, set up by municipalities or to a dealer providing this service. The separate disposal of the product avoids possible negative consequences for the environment and human health, deriving from inappropriate disposal and it allows you to recover the materials it is made in order to obtain significant savings in energy and resources. To emphasize the need to dispose of medical equipment separately, the product is marked with

the crossed out wheeled bin.

TREATMENT OF USED BATTERIES - (Directive 2006/66/EC) 

This symbol on the product indicates that the batteries should not be treated as household waste. Making sure the batteries are properly disposed of will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate disposal. The recycling of materials will help to protect natural resources. Give the batteries to the collection points for recycling. For more detailed information about the disposal of batteries or product, please contact your municipality, waste disposal service or the shop where you purchased the unit.

TECHNICAL FEATURES

LTK600

Dimensions	330 x 330 x 25 mm
Weight	1 Kg
Max. loading weight	180 Kg
Min. resolution	100 g
Display	Digital LCD 4"
Supply	4 batteries 1,5 V type AAA

Special features

- Switching on by foot contact
- Low battery and overload indicator
- Auto reset and power off
- Units to be chosen among Kg/Lb/St

LTK610

Dimensions	300 x 300 x 30 mm
Weight	1,5 Kg
Max. loading weight	150 Kg
Min. resolution	100 g
Display	Digital LCD 2,5"
Supply	1 lithium battery 3V model CR2032

Special features

- High-precision sensor
- Units to be chosen among Kg/Lb/St
- Low battery and overload indicator

LTK615

Dimensions	Diameter 330 mm x 20 mm high
Weight	1,5 Kg
Max. loading weight	150 Kg
Min. resolution	100 g
Display	Digital LCD 3"
Supply	1 lithium battery 3V model CR2032

Special features

- High-precision sensor
- Switching on by foot contact
- Units to be chosen among Kg/Lb/St
- Low battery and overload indicator
- Auto reset and power off

WARRANTY

Moretti products are guaranteed from material or manufacturing faults for 2 years from the purchasing date, except possible exclusion or restriction as follows. The warranty shall not be apply in the event damages caused by improper use, abuse or alteration, and for the warranty will not be valid if the instructions for use are not strictly adhered to. The correct intended use is specified in this manual.

Moretti is not responsible of consequent damages, personal injuries or whatever caused by or in relation with mistaken installation or with device improper use.

Moretti does not guarantee for damages caused by the following condition: natural disaster, not authorized maintenance or repairs, faults caused by problems on electricity supply (when necessary), use of spare parts not covered by Moretti, improper use, alteration not authorized, shipment damages (different from original Moretti shipment), or in case of insufficient maintenance as indicated in the manual.

The warranty doesn't cover components subject to wear and tear during the correct use of the device.

RECHARGEABLE BATTERY WARRANTY (IF INCLUDED)

The original and rechargeable batteries are under a 90 days warranty referring to performances and a 6 month warranty referring to production defects or according to the law requests.

If charged batteries are not used for more than 3 months, the warranty lapses.

If law batteries are not used for more than 3 days, the warranty lapses.

REPAIRING

Warranty repair

If a Moretti item presents material or manufacturing faults during the warranty period, Moretti will estimate with customer if the fault can be cover from warranty. Moretti, at its unquestionable discretion, could repair or replace the item, with a Moretti dealer or with own base. Labor cost can be in charge to Moretti if the repair is cover from warranty. A repair or a replacement doesn't extend the warranty.

Repair not cover from warranty

A product without warranty can be given back after Moretti authorization. The labor and shipping costs for a bath lift out of warranty are to be payed by the customer or by the dealer. The repairs are guaranteed for 6 months as from the good received.

Non-defective device

The customer will be informed if, after the device return and examination, Moretti declares that the device is not faulty. In this case the good will be sent back to customer, the shipping charge will be to be payed by the customer.

Spare parts

The original Moretti spare parts are guaranteed for 6 months as from the receiving date.

EXEMPT CLAUSES

Moretti doesn't offer any other declarations, warranty or conditions, explicit or implicit, including possible declarations, warranties or condition of merchantability, fitness for a specific purpose, non infringement and non interference, all but what expressly specified in this warranty. Moretti doesn't guarantee the nonstop and faultless usage. The duration of possible implicit warranties which can be fixed by according to the law is limited by the warranty period, in the limits of law. Some states or countries don't allow to limit the implicit warranty or the exclusion or the limitation for accidental damages. In such countries, some of those exclusions or limitation may not be applied to the user. The present warranty may be modified without prior notification.

WARRANTY CERTIFICATE

Product _____

Purchasing date _____

Authorized dealer _____

Street _____ **Place** _____

Sold to _____

Street _____ **Place** _____

Via Bruxelles, 3 - Meleto 52022 Cavriglia (Arezzo) Tel. +39 055 96 21 11

www.morettispa.com email: info@morettispa.com

MADE IN P.R.C.

MORETTI S.p.A.
Via Bruxelles 3 - Meleto
52022 Cavriglia (Arezzo)

Tel. +39 055 96 21 11
Fax +39 055 96 21 200

www.morettispa.com
info@morettispa.com